

## CELLULAR MATRIX™



“THE PAIN AND MOBILITY  
IMPROVEMENTS OF HA BIOPOLYMER  
COMBINED WITH THE CELL  
STIMULATION PROPERTIES AND  
ANTI-INFLAMMATORY ACTION  
OF A-PRP GROWTH FACTORS”

[WWW.CELLULARMATRIX.INFO](http://WWW.CELLULARMATRIX.INFO)




A-PRP & HA ALL IN ONE

regenlab

# CELLULAR MATRIX™

“THE PAIN AND MOBILITY IMPROVEMENTS OF HA BIOPOLYMER COMBINED WITH THE CELL STIMULATION PROPERTIES AND ANTI-INFLAMMATORY ACTION OF A-PRP GROWTH FACTORS”

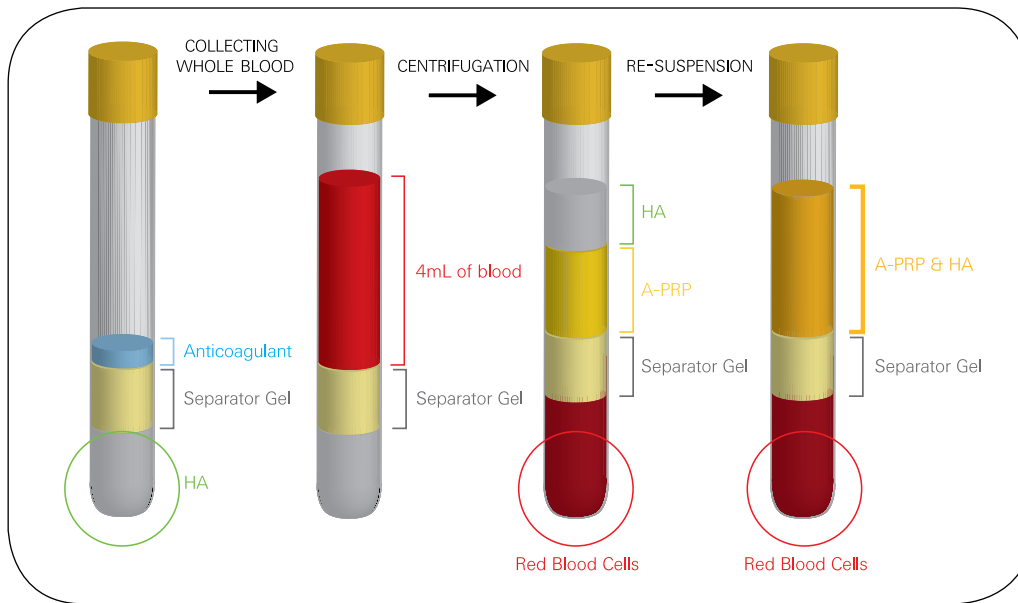

## ***Synergistic action in the treatment of osteoarthritis***


### HYALURONIC ACID

- Major component of synovial fluid contributing to joint homeostasis
- Molecular weight and concentration of HA decrease in OA
- Viscosupplementation provides pain relief and improved function with 1-3 i.a. injections of HA
- 25 years of clinical experience show benefits of HA last 6-12 months with 75% response rate
- Network of HA chains constitutes ideal cell-friendly matrix (Gobbi 2009)

### AUTOLOGOUS PLATELET-RICH PLASMA

- Autologous reservoir of growth factors from patient's own blood
- Proven role in the healing of soft and hard tissues (Marx 2004)
- Growing body of literature in the treatment of OA with 1-3 i.a. injections in studies with 6-12 month follow-up
- Mechanism of action comprises anti-inflammatory activity and activation of cell-signalling cascades (Barry 2001, Sevens 2004)
- Key role in new matrix synthesis for tissue regeneration

# ONE STEP PROCEDURE ; CLOSED SYSTEM

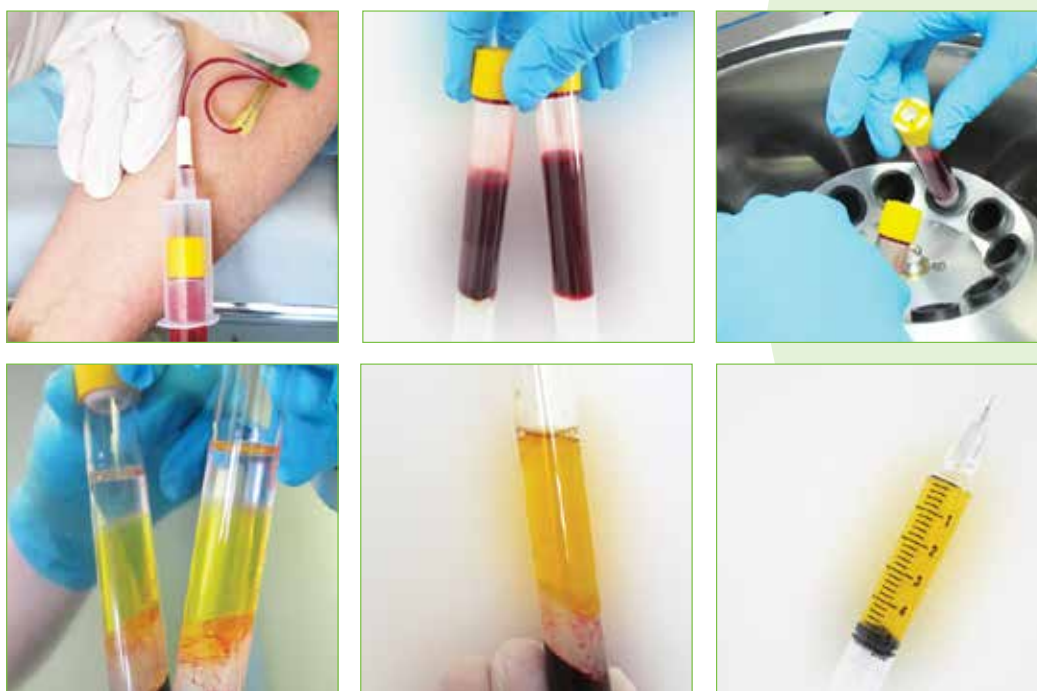

## CELLULAR MATRIX™ PRODUCT CHARACTERISTICS

- Contains 2ml of natural, non-crosslinked, fermented HA (1,550 kDa) at a concentration of 20 mg/ml (40 mg total) in addition to the thixotropic cell-separation gel
- Provides 4ml of the A PRP-HA combination in a closed system with a single gesture from 4 ml of autologous blood
- Immediate formation of a cell-friendly HA network in which platelets are dispersed
- Significantly higher viscosity than A-PRP
- Potential for increased residence time and growth factor release in joint
- Positive impact on cell migration and differentiation
- Treatment regimen of up to 3 i.a. injections
- 3 standard Regen pharmaceutical-grade glass tubes packaged in double blister, sealed with 2FS Tyvek® in accordance with ISO 11607



INTENDED USE OF THE DEVICE:  
**Intra-articular injections for symptomatic treatment  
of articular pain and mobility improvement.**

PATENTED BY REGEN LAB SA

International Publication Number

**WO2011/110948**



CELLULAR MATRIX™



ACCESSORY KIT



regenlab



REGEN LAB SA

En Budron B2

CH - 1052 – Le Mont-Sur-Lausanne (CH)

Tel. +41(0)21 864 01.11

Fax. +41(0)21 864 01.10

info@regenlab.com

www.regenlab.com